

SafeMask® Premier Elite

Procedure Earloop Face Masks



ASTM LEVEL 3

2047

2044

2043

2042

2046


50 masks/box; 10 boxes/case Also available in ASTM levels 1 & 2.





SafeMask® NEW Curve Procedure Earloop Face Masks **Contoured fit ensures** breathability and more complete coverage. **ASTM LEVEL 3** 220114 220115 25 masks/box; 8 boxes/case



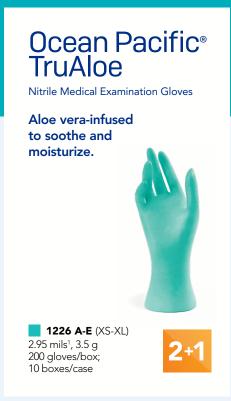



SafeMask® SofSkin

Procedure Earloop Face Masks

Welcome to the SafeMask SofSkin collection.



SafeMask® SofSkin ClearFit

Skin-friendly material + a chin piece + integrated fog-free strip to prevent fogging of eyewear and loupes.

ASTM LEVEL 3

220014 White

50 masks/box; 10 boxes/case

SafeMask® SofSkin Fit

Skin-friendly material + a chin piece providing an even more personalized fit that conforms to each individual face shape.

ASTM LEVEL 3

219914 White

50 masks/box; 10 boxes/case

SafeMask® SofSkin

Skin-friendly, ultra-soft, and dye-free material made especially for sensitive skin.

ASTM LEVEL 3

2086 White

2088 Pink

2087 Blue

2089 Lavender

50 masks/box; 10 boxes/case Also available in ASTM level 1.



SafeMask® Classics

Procedure Earloop Face Masks

Reliable respiratory protection at a suprisingly low price.



205316

50 masks/box; 10 boxes/case Also available in ASTM levels 1 & 2.

Ocean Pacific® NeoNatural

Polychloroprene Medical Examination Gloves

Comfort and strength.



1228 A-E (XS-XL)

4 mils¹, 6 g 100 gloves/box; 10 boxes/case

Ocean Pacific® Elements

Biodegradable Nitrile Medical **Examination Gloves**

81% biodegradable in less than 1.5 years.*



OP(XS-XL)-BIO

3 mils¹, 3.5 g

200 gloves/box; 10 boxes/case



Duraflor Ultra®

White 5 % Sodium Fluoride Varnish

Fast-acting, great tasting, dries clear.

Effective treatment of hypersensitivity. Dentist and patient preferred ready-to-use and easy-to-apply, nonclumping formula. Does not contain gluten, peanuts, tree nuts, soy, sesame, corn, shellfish, egg or dairy products.

1016-ST30 Strawberry **1016-CH30** Cherry

1016-M30 Mint

1016-BG30

Bubble gum

0.4 ml x 30 units/box Promotion also applicable on 200 units/box



HemoStat™

Topical Hemostatic Dressing

Control bleeding stat. No suturing required.



SafeMatrix™

Single-Use Matrix Band





20 units/box; 200 boxes/case



SafeSeal[™] Quattro

Self-Sealing Sterilization Pouches

Cover all your bases

Multivariable process indicators confirm time, temperature and pressure.

 88000-4
 2 ¾ x 9"
 88025-4
 5 ¼ x 10"

 88005-4
 3 ½ x 5 ¼"
 88030-4
 7 ½ x 13"

 88010-4
 3 ½ x 9"
 88035-4
 10 x 14"

 88015-4
 2 ¼ x 4"
 88040-4
 12 x 17"

 88020-4
 5 ¼ x 6 ½"
 88045-4
 4 ¼ x 11"

200 pouches/box; 6 boxes/case

88010-4B 3 ½ x 9"

500 pouches/box; 6 boxes/case





SafeWear™ High Performance Lab Coat

Premium Protective Apparel

Protection and style featuring knitted cuffs and collar.



- **8107 (A-F)** Deep Blue
- 8108 (A-D) Pretty Pink
- 8109 (A-D) Plum Purple
- 8110 (A-D) White Frost
- **8112 (A-D)** Soft Blue
- 8117 (A-D) Tropical Teal
- **8119 (A-D)** Jet Black

12 units/bag

Also available in hip-length jacket.

Hedy™ Dental Dam

Polyisoprene

Elongation up to 1230%. Premium, tear-resistant, non-latex, powder-free.

320DNB-5M 5 x 5" - Medium

320DNB-5H 5 x 5" - Heavy

320DNB-6M 6 x 6" - Medium

320DNB-6H 6 x 6" - Heavy

30 units/box Also available in latex.



SafeBasics® Headrest Covers

Tissue/Poly

Reduce the risk of crosscontamination. Fluid repellent and breathable. Protect chair materials.

3015 Blue

3016 Dusty Rose

3020 Lavender

3014 White

 10×10 " | 500 units/case Also available in 10×13 ".





FOR THOSE WHO ENSURE HEALTH AND CREATE BEAUTIFUL SMILES

Our proven, high quality, single-use infection prevention and control solutions for PPE, sterilization, office essentials, and preventive care are thoughtfully designed for simplicity and comfort without compromise to safety and clinical efficacy.

Making each day in your practice a little easier.

Founded in 1988 in response to the critical demand for high-quality medical gloves during the global HIV crisis, Medicom has consistently been at the forefront of infection prevention and control innovation.

Drawing on over 35 years of expertise, we have continuously evolved our mask technology and designs to better meet the wide array of individual needs and preferences.

Explore our extensive mask portfolio today and experience the difference that quality, innovation, and fit can make in ensuring your safety and comfort.

CHOOSING THE RIGHT MASK

Protection and comfort

The ASTM standard F2100-23 classifies masks into three levels based on their bacterial filtration efficiency (BFE), sub-micron particulate filtration efficiency (PFE), fluid resistance and breathability (differential pressure).

	Low Barrier ASTM Level 1	Moderate Barrier ASTM Level 2	High Barrier ASTM Level 3
BFE, at 3 microns	≥ 95%	≥ 98%	≥ 98%
PFE, at 0.075 micron	≥80%	≥ 85%	≥ 85%
Fluid Resistance	≥80 mm Hg	≥120 mm Hg	≥160 mm Hg
Differential pressure	< 5.0 mm H ₂ O/cm ²	< 6.0 mm H ₂ O/cm ²	< 6.0 mm H ₂ O/cm ²
Flammability	Class 1	Class 1	Class 1

Determine the level of barrier protection you need, taking into account factors such as procedure duration, amount of fluid or aerosol involved, and standard precautionary measures.

ASTM LEVEL 1*

- Low-risk exams
- Operatory cleaning/maintenance
- Orthodontics
- X-rays, impressions
- Lab trimming, finishing, and polishing
- Fluoride application

ASTM LEVEL 2*

- Limited oral surgery
- Endodontics
- Prophylaxis
- Bleaching/whitening treatments
- Restoratives/composites
- Sealants
- Scaling and root planing

ASTM LEVEL 3*

- Complex oral surgery
- Crown preparation
- Periodontal procedures
- Implants
- Using ultrasonic scalers
- Using high-speed hand piece
- Laser based applications

Form and function

Medicom offers a variety of general purpose, procedure and surgical masks in all three ASTM levels, and NIOSH-certified N95 surgical respirators. Our masks are not made with natural rubber latex and will not lint or tear during use.



Surgical Tie-on Face Mask

For dental/clinical procedures.



Surgical N95 Respirator



For treating patients with known or suspected aerosol transmissible diseases.

ANATOMY OF A FACE MASK

Filter layer Earloops/ties o Adjustable nose piece **Dye-free** inner layer Outer layer ·······o **Adjustable** chin piece Construction





Shingle pleats



Omega pleats (larger air pocket)



3D fit



Window

Specialty considerations

3D Fit



Greater coverage Window



Speechreading enabled

Larger Air **Pocket**



Better breathability **Innovative Filter**



Lightweight

Tangle-Free Dispensing



Workflow efficiency

Skin Friendly



Designed for sensitive skin

Moisture Wicking



Prevents moisture build-up

Ultra-Soft Inner Layer



Added comfort

Fog-Free Strip



Prevents fogging

Adjustable Chin Strip



Personalized fit

MASKS AS INDIVIDUAL AS YOU



SafeMask® Curve



ASTM Level 3	
220114	
220115	

25/box, 8 boxes/case, 200/case





ASTM Level 3

221415 221492

221414

SafeMask® Artisan

ASTM Level 2

221315



Summer 2025





SafeMask® Vision





ASTM Level 3	
221514	
50/hox 10 hoxes/case 500/case	







ASTM Level 3

2086

2087

2088 2089



SafeMask® SofSkin Fit





SafeMask® SofSkin



ASTM Level 1

2080

2081

2082

2083



Le	vel	3

ASTM Level 3	
219914	
50/box, 10 boxes/case, 500/case	







ASTM Level 3	
220014	
E0/hay 10 hayaa/aaaa E00/aaaa	

50/box, 10 boxes/case, 500/case











SafeMask® TailorMade





50/box, 10 boxes/case, 500/case







ASTM Level 1	ASTM Level 2	ASTM Level 3
200311		200511
200315	200415	200515
200316	200416	200516
		200514

50/box, 10 boxes/case, 500/case

SafeMask® Master Series



ASTM Level 1	ASTM Level 2	ASTM Level 3
2050		
2051	2056	2058
2052		
2054		2059
	2055	2057

50/box, 10 boxes/case, 500/case













SafeMask® Premier (Family)



SafeMask® Premier Elite ProShield

Midnight

ASTM Level 1	ASTM Level 2	ASTM Level 3	
2010		2047	
2011		2044	
2015	2040	2042	
2016		2046	
2017			
2018		2043	

2010		2047	
2011		2044	
2015	2040	2042	
2016		2046	
2017			
2018		2043	

ASTM Level 3 204215 204216 2025 25/box, 4 boxes/case, 100/case

ASTM Level 1	ASTM Level 2	ASTM Level 3	
213792	213892	213992	
50/box, 10 boxes/case, 500/case			

50/box, 10 boxes/case, 500/case











SafeMask® Classics

ASTM Level 1	ASTM Level 2	ASTM Level 3
205115	205215	205315
205116	205216	205316
205117	205217	205317

50/box, 10 boxes/case, 500/case

205214

AssureMask® Precision





ASTM Level 1	ASTM Level 2	ASTM Level 3
206815	207015	207115
206818	207018	207118
206814	207014	207114

50/box, 10 boxes/case, 500/case

SafeMask® Tie-On

ASTM Level 1	
2000	

50/box, 6 boxes/case, 300/case



SafeMask® Junior



205314





SafeMask® Cone

ASTM Level 1	
2035	

50/box, 10 boxes/case, 500/case

2020

50/box, 12 boxes/case, 600/case




FIND YOUR MATCH

Procedure Masks	ASTM protection level	Construction*	Innovative filter	Skin friendly	Moisture wicking	Fog-free strip	Adjustable chin strip	Ultra-soft inner layer	Colors
SafeMask® Curve	L3	3D fit							• 0
SafeMask® Artisan	L2, L3		•						0 • •
SafeMask® Vision	L3	Window	•						0
SafeMask® SofSkin	L1, L3			•	•				0 • • •
SafeMask® SofSkin Fit	L3			•	•		•		0
SafeMask® SofSkin ClearFit	L3			•	•	•	•		0
SafeMask® TailorMade	L1, L3						•		•
SafeMask® FreeFlow	L1, L2, L3	Omega pleats				•			0 • • •
SafeMask® Master Series	L1, L2, L3							•	•••
Midnight	L1, L2, L3								•
SafeMask® Premier (Family)	L1, L2, L3								• 0 •
SafeMask® Premier Elite ProShield (with visor)	L3								• • •
SafeMask [®] Classics	L1, L2, L3								0 • • •
SafeMask® Junior (pediatric)	L1								•

Other Masks	ASTM protection level	Construction*	Innovative filter	Skin friendly	Moisture wicking	Fog-free strip	Adjustable chin strip	Ultra-soft inner layer	Colors
AssureMask® Precision	L1, L2, L3					•			••0
SafeMask® Tie-On	L1								•
SafeMask® Cone	N/A	Cone							•

^{*}Shingle pleats, unless indicated otherwise.

SURGICAL N95 RESPIRATORS

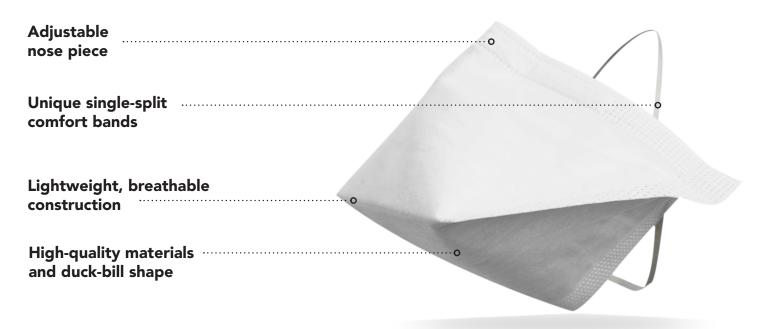

SafeMask Architect Pro delivers the highest standards protection without compromise to comfort.

Test Method	Results			
Biocompatibility: cytotoxity, skin sensitization and skin irritation ISO 10993-5, ISO 10993-10 and ISO 10993-23	Complies with the standard			
Particulate filter efficiency (PFE) 42 CFR 84	≥ 95 %			
Airflow resistance 42 CFR 84	Inhalation < 35 mm H ₂ O Exhalation < 25 mm H ₂ O			
Fluid resistance (resistance to penetration by synthetic blood) ASTM F1862	Level 2: 120 mm Hg			
Flammability 16 CFR 1610	Class 1 - Complies			
Fit CAN/CSA Z94.4	A combination of all sizes fits > 90% of users			

Danilla

Regulations require respirator users to be fit-tested to confirm that the respirator forms a tight seal around the face before using it in the workplace.



Scan to learn more about fit testing and/or to access step-by-step videos, sizing guides, donning and doffing instructions and more.

SafeMask® Architect Pro™ NIOSH-certified respirator that filters out 95% of particulate aerosols free of oil 203214 Small 203314 Medium 203414 Large







Dental USA medicom.com | 1-800-361-2862

Rev.02 March 2025 ©2025 AMD Medicom Inc. All rights reserved.



Reliable protection without compromise to comfort and style

Reduce the risk of cross-contamination and eliminate the need for laundering with our portfolio of high-quality, single-use PPE options.

SafeBasics™

SafeBasics protective apparel provides a basic level of protection at an attractive price.

- Combines protection and comfort
- Versatile design and material options
- Cost effective



Single-use gown for procedures requiring protection from spray and splatter. Fluid-repellent design provides added protection against cross-contamination.



Lab Gowns

Front snap closure, knitted cuffs and collar for comfortable, easy on-and-off fit.



Fluid Impervious Gowns

Polyethylene coating provides fluid resistance for procedures involving the potential for crosscontamination.



Item	Description	Unit/ Pkg	Pkg/ Case
D-8009	Isolation Gown Reg - Yellow	10	5
D-8010	Isolation Gown XL - Yellow	10	5
8023	Isolation Gown Reg - Blue	10	5

Item	Description	Unit/ Pkg	Pkg/ Case
8016	Lab Gown S - Sky Blue	10	5
8017	Lab Gown M - Sky Blue	10	5
8018	Lab Gown L - Sky Blue	10	5
8019	Lab Gown XL - Sky Blue	10 e	5

Item	Description	Unit/ Pkg	Pkg/ Case
8022	Fluid Impervious Gown Reg - White	10	5
8024	Fluid Impervious Gown XL - White	10	5



Bouffant Caps

Hair and head fluid protection with breathable, low-linting design.

Beard Covers

Face protection with an elastic edge and headband loop for optimal fit and comfort.

Item	D	escription	Unit/ Pkg	Pkg/ Case
8000		Bouffant Caps 21" - White	100	5
8001		Bouffant Caps 21" - Blue	100	5
8003		Bouffant Caps 24" - White	100	5
8004		Bouffant Caps 24" - Blue	100	5
8035		Beard Covers - White	100	10



Fluid repellent footwear protection.



Item	D	escription	Unit/ Pkg	Pkg/ Case
8006		Shoe Cover Non-Conductive Non-Skid - XL - Blue	100	3
8007		Shoe Cover Non-Conductive Non-Skid - Reg - Blue	100	3
8008		Shoe Cover Non-Conductive Reg - Blue	100	3
8036		Shoe Cover Non-Conductive Non-Skid - XL - White	100	3

SafeWear™

Premium Protective Apparel

- Lightweight, soft, supple fabrics keep wearers cool and comfortable.
- Secure front snaps.
- Latex free.
- Made from AAMI PB70 Level 2 tested material.
- Complies with recommended OSHA best practices.
- Available in a wide range of sizes, from XS 3XL.





Form-Fit Isolation Gown™
Premium isolation gown

Item	D	escription	Unit/ Pkg	Pkg/ Case
8113¹		Sunny Yellow	12	5
8114¹		Bright Blue	12	5



Hipster Jacket[™] Premium hip length jacket

Item	D	escription	Unit/ Pkg	Pkg/ Case
8101 ²		Deep Blue	12	5
8102*		Pretty Pink	12	5
8103 ²		Plum Purple	12	5
8104		Soft Blue	12	5
8105		White Frost	12	5
8115*		Tropical Teal	12	5
8116*		Poppy Pink	12	5
8118*		Jet Black	12	5



High Performance Lab Coat™ Premium lab coat

Item	D	escription	Unit/ Pkg	Pkg/ Case
8107		Deep Blue	12	5
8108*		Pretty Pink	12	5
8109		Plum Purple	12	5
8110		White Frost	12	5
8112		Soft Blue	12	5
8117*		Tropical Teal	12	5
8119*		Jet Black	12	5

¹Available only in regular size.

Size:						
XS: Extra Small	A: Small	B: Medium	C: Large	D: Extra Large	E: 2X Large	F: 3X Large
² Available only in Hipster Jacket, Deep blue & Plum Purple.	Hipster Jackets a 2X and 3X Large	Available in US only.				

^{*}Colour available in US only.

Protective Apparel Best Practices

DONNING

- Perform hand hygiene before removing item from package.
- Select appropriate garment for intended use (length, coverage, fluid resistance, etc.)
- Fasten correctly.

DOFFING

- Remove apparel immediately after use and dispose in a designated receptacle. Do not re-use.
- Perform proper hand hygiene after removing apparel to prevent cross-contamination.

KEEP IN MIND:

- Protective apparel should not be worn in restrooms, personal offices or break rooms.
- Change protective apparel between patients and as soon as possible when penetrated by blood or other potentially infectious fluids.
- Turn garment inside out when removing to keep infectious material inside.



medicom.com | 1-800-361-2862

